



Patient Last Name: \_\_\_\_\_ Patient First Name: \_\_\_\_\_  
 Fitter Last Name: \_\_\_\_\_ Fitter First Name: \_\_\_\_\_  
 Fitter Title: \_\_\_\_\_ (example: PT/OT/PTA)  
 Date: \_\_\_\_\_

## MEDI NON-CUSTOM CIRCULAR KNIT ARM SLEEVE MEASUREMENT FORM AND SIZING CHART

### Measure in centimeters

Left  Right

**Length Measurement from axilla to wrist (c-g) = \_\_\_\_\_**



The length measurement is taken from **point 'c'** along the inside of the arm to the inseam, **point 'g'**, of the arm; measure with patient sitting and arm extended straight on table.

### Measurement Points

Circumference **point 'g'** is taken around the upper arm at the axilla. For ease in determining this point, place a folded piece of paper up into the arm pit area. Point 'g' on the arm will be even with the top of the paper. **Measure with moderate tension.**

Circumference **point 'e'** is found in the middle of the elbow. **Measure without tension.**

Circumference **point 'c'** is the point of greatest compression and therefore the most important point. This point is found distally on the arm where the wrist meets the hand. One should find the smallest circumference at this point. **Measure without tension.**

**Circumference 'g'** = \_\_\_\_\_ **c-g with silicone border**

**Circumference 'e'** = \_\_\_\_\_ **c-g with no silicone**


**Circumference 'c'** = \_\_\_\_\_

**Circumference 'c1'** = \_\_\_\_\_

Item	Class I	Class II	Colors Available
Harmony	20-30 mmHg <input type="checkbox"/>	30-40 mmHg <input type="checkbox"/>	Sand <input type="checkbox"/> Caramel <input type="checkbox"/> Black <input type="checkbox"/>

### SIZE CHART

## lymphology measuring charts

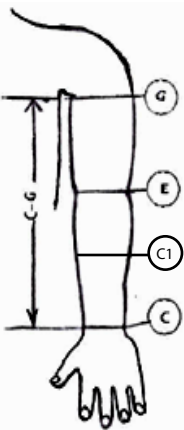

mediven® harmony sleeves, sleeve with gauntlet, gauntlet and gloves

circumference (cm)	1	2	3	4	5	6	7	8
"G" extra-wide	28-31	30-34	32.5-36.5	34.5-39.5	36.5-41.5	39-44	41-46.5	44.5-50.5
"G" standard	25-28	27.5-31	29.5-33.5	31.5-36	33-37.5	35.5-40.5	37.5-42.5	40-45.5
"E" extra-wide	24-26.5	26-29.5	28-31.5	29.5-34	31-35.5	33-37.5	35-40	38-43
"E" standard	23-25	24-27	25.5-29	27-31	28.5-32.5	30-34.5	32-36.5	34.5-39.5
"C1"	15-16.5	16-18.5	17.5-20	19-22	20.5-24	22.5-26	24.5-28.5	27-31
"C"	13.5-14.5	14-16	15.5-18	17-20	18.5-21.5	20-23.5	22-25.5	24.5-28.5
"A"	15-18	15.5-19.5	16.5-20.5	18.5-23	19.5-24	21.5-26	23.5-28	25.5-30

gloves & gauntlets available in sizes 2-7

armsleeves available in sizes 1-8