



Patient Last Name: _____ Patient First Name: _____
 Fitter Last Name: _____ Fitter First Name: _____
 Fitter Title: _____ (example: PT/OT/PTA)
 Date: _____



JUXTAFIT LOWER LEG GARMENT MEASUREMENT FORM

SIZING CHART


Patient should be standing when taking measurements. Limbs should be at their smallest; limb size should be maintained from time of measurements to ensure a proper fit of garment. When on borderline between sizes, select a smaller size if significant edema is present and reduction is anticipated; if edema is minimal, select larger size. After measuring the smallest dimension of the ankle and the largest dimension of the calf, if any of the measurements do not fall within one of these size ranges, call (800) 247-2243. Length measurement is for Juxta-Fit™ only, actual legging lengths are 28cm (Short) and 36cm (Long); Ready-Fit™ and Classic-Flex™ have one adjustable length. **Note 1in = 2.54cm**

MEASUREMENTS		
	LEFT	RIGHT
B		
B₁		
C		
H		
B-D		

JUXTA-FIT™ PREMIUM & ESSENTIALS			
SIZE	CIRCUMFERENCE MEASUREMENT		
	B	B ₁	C
SMALL	15-23cm	20-28cm	30-38cm
MEDIUM	20-28cm	25-33cm	36-43cm
MEDIUM FULL CALF	20-28cm	33-40cm	43-50cm
LARGE	25-33cm	30-38cm	41-48cm
LARGE FULL CALF	25-33cm	38-46cm	48-56cm
XL	30-38cm	33-41cm	41-51cm
XL FULL CALF	30-38cm	41-51cm	51-61cm
XXL*	33-41cm	43-53cm	53-64cm

AFW MEASUREMENT	
SIZE	H
SMALL	19-24cm
MEDIUM	24-29cm
LARGE	29-34cm

LENGTH MEASUREMENT	
SIZE	B-D
SHORT	under 36cm
LONG	over 36cm

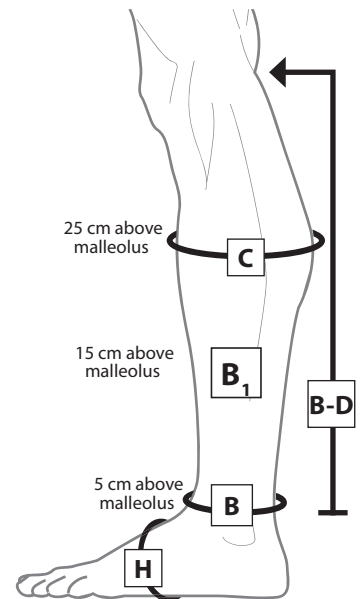

SIZE OF JUXTAFIT LEG NEEDED: _____
 SIZE OF JUXTAFIT AFW NEEDED: _____
 QTY OF JUXTAFIT LEG NEEDED: _____
 QTY OF JUXTAFIT AFW NEEDED: _____

* If measurements fall within XXL size range, please consider ordering a custom Juxta-Fit™ for optimal fit.
 ‡ EZ Single-Band Ankle-Foot Wrap™ included by default
 § Juxta-Lite™ Ankle-Foot Wrap is available only in sizes Small and Medium
 † Comfort Compression Anklets™ included by default

All products contain the CircAid Built-In-Pressure (BPS) guide card. Product lengths may differ from body measurements.
 Size SMALL JuxtaFit receive STANDARD Compression Ankle. All other sizes receive LARGE Compression Ankle.
 Size SMALL, MEDIUM, MEDIUM-FULL CALF and LARGE JuxtaFit receive SMALL Cover Up. All other sizes receive LARGE Cover Up.
 Size SMALL, MEDIUM and MEDIUM-FULL CALF JuxtaFit receive STANDARD Pac Band. All other sizes receive LARGE Pac Band.