



Patient Last Name: _____ Patient First Name: _____
 Fitter Last Name: _____ Fitter First Name: _____
 Fitter Title: _____ (example: PT/OT/PTA)
 Date: _____

CIRCAID CUSTOM ARMSLEEVE
 SIZING CHART
Making Compression Easier.

The Juxta-Fit™ Armsleeve is constructed from the wrist to the axilla measurement. A lateral rise is incorporated at the top of the garment in addition to a 10° bend is at the elbow. Measure patient's arm length with the arm bent at 35°. Also ensure that the arm is at its smallest and maintain arm size at time of measuring to assure a proper fit of the garment(s). **NOTE: 1in = 2.54cm**

ARM AND HAND SELECTION	
<input type="checkbox"/> RIGHT	<input type="checkbox"/> LEFT

COVERUP COLOR
<input type="checkbox"/> BLACK
<input type="checkbox"/> BEIGE
<input type="checkbox"/> PINK

JUXTA-FIT™ GAUNTLET SIZE OPTIONS					
PALM CIRCUMFERENCE	14-17cm	17-20cm	20-23cm	23-26cm	26-29cm
SIZE	<input type="checkbox"/> XSMALL	<input type="checkbox"/> SMALL	<input type="checkbox"/> MEDIUM	<input type="checkbox"/> LARGE	<input type="checkbox"/> XL

ARM LENGTH	
c-e	c-g
wrist to elbow crease	wrist to axilla (or where the garment should end)

HAND OPTIONS
<input type="checkbox"/> GAUNTLET
<input type="checkbox"/> GLOVE
<input type="checkbox"/> DORSUM STRAP
<input type="checkbox"/> NONE

JUXTA-FIT™ GLOVE MEASUREMENTS										
LENGTH	1	2	3	4	5	6	7	8	9	10
	1x	2x	3x	4x	5x	a	b	c		
CIRCUMFERENCE	1w	2w	3w	4w	5w					

ARM CIRCUMFERENCE MEASUREMENTS (Palm; Wrist; 5cm increments from the wrist to the axilla)	MEASURING POINTS
	Palm
	Wrist
	5cm
	10cm
	15cm
	20cm
	25cm
	30cm
	35cm
	40cm
	45cm
	50cm
	55cm


MEASUREMENTS 1w-5w: the circumference where the glove should end; commonly, the Distal Interphalangeal joint (DIP).

MEASUREMENTS 1x-5x: the circumference of each finger at the Proximal Interphalangeal joint (PIP) or if the base of the finger is larger than the PIP please use the larger measurement.

MEASUREMENT a: circumference of palm at base of pinky finger

MEASUREMENT b: circumference of palm at base of thumb webbing

MEASUREMENT c: circumference of the wrist

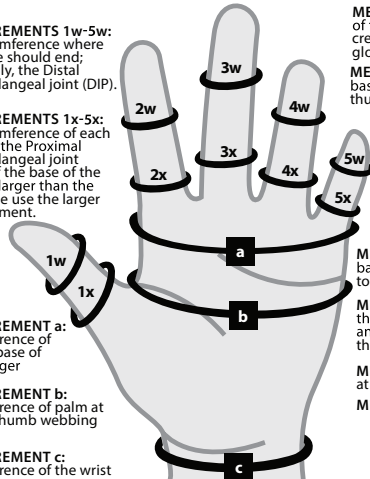

MEASUREMENTS 1-5: length of finger from the center of the crease at base to where the glove should end

MEASUREMENT 6: length from base of index finger to the base of thumb webbing

MEASUREMENT 7: length from base of index finger (same as 6) to base of palm

MEASUREMENT 8: length from the crease between the ring finger and pinky finger to the base of the palm

MEASUREMENT 9: width of palm at base of thumb webbing

MEASUREMENT 10: width of palm at base of pinky finger

