



Patient Last Name: _____ Patient First Name: _____
 Fitter Last Name: _____ Fitter First Name: _____
 Fitter Title: _____ (example PT/OT/PTA)
 Date: _____

Measure in Centimeters

Left Right

JUZO NON-CUSTOM CIRCULAR KNIT LEG STOCKING MEASUREMENT FORM AND SIZING CHART

Thigh Length _____

Calf Length _____

A-g REGULAR LENGTH IS 72-83 CM
 A-g SHORT LENGTH 63-72 CM
 A-g PETITE LENGTH 55-63 CM

A-d REGULAR LENGTH IS 40-46 CM
 A-d SHORT LENGTH IS 34-40 CM
 A-d PETITE LENGTH IS <33 CM

For thigh length, the length measurement is taken from the floor (at the inside of heel) along the inside of the leg to the inseam/groin, point 'g;' this length should end at the gluteal fold.

For calf length, the length measurement is taken from the floor (at the inside of heel) along the inside of the leg to approximately 3 1/2 centimeters below the crease of the knee, point 'd'

Silicone borders are not usually recommended for calf length stockings, however, they shaped legs.

For calf and thigh length stockings, measure with the patient standing.

Measurement Points

Circumference **point 'g'** is found at the inseam/groin

Circumference **point 'c'** is the largest circumference of the calf

Circumference **point 'b'** is the point of greatest compression and therefore the most important point. This point is the smallest circumference of the ankle found closest to and above the ankle bone

Circumference **'g'** = _____

Circumference **'c'** = _____

Circumference **'b'** = _____

Open Toe Silicone

Closed Toe No Silicone


Measurement Point for the Body Part (KM) or Hip Attachment

*The following measurement will be necessary if you are ordering a "custom-made" thigh length stocking to be sewn to a body part.

*It is common for a Class 1 Body Part (KM) to be sewn to a Class 2 or 3 leg stocking.

Hip circumference = _____

KM Class 1 (18-21mmHg) KM Class 2 (23-32mmHg)

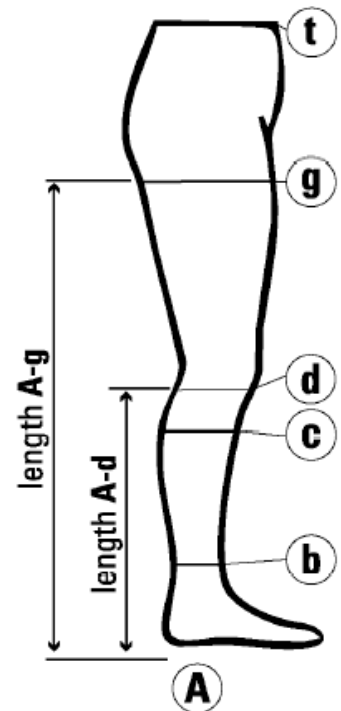

Material <small>*Dynamic is new name for Varin</small>	Class 1 20-30 mmHg	Class 2 30-40 mmHg	Class 3 40-50 mmHg	Colors Available
Dynamic (calf length) AD	3511 <input type="checkbox"/>	3512 <input type="checkbox"/>	3513 <input type="checkbox"/>	Beige <input type="checkbox"/> Chestnut <input type="checkbox"/> Black <input type="checkbox"/>
Dynamic (thigh length) AG	3511 <input type="checkbox"/>	3512 <input type="checkbox"/>	3513 <input type="checkbox"/>	Beige <input type="checkbox"/> Chestnut <input type="checkbox"/> Black <input type="checkbox"/>
Dynamic Full Knit (thigh length) AG		3512L <input type="checkbox"/>	3513L <input type="checkbox"/>	Beige <input type="checkbox"/> Chestnut <input type="checkbox"/> Black <input type="checkbox"/>
Dynamic Full Knit (panty hose)			3513ATLY <input type="checkbox"/>	
Soft 2000 (calf length) AD	2001 <input type="checkbox"/>	2002 <input type="checkbox"/>		Beige <input type="checkbox"/> Chestnut <input type="checkbox"/> Black <input type="checkbox"/> _____ <input type="checkbox"/> * Seasonal colors available
Soft 2000 (thigh length) AG	2001 <input type="checkbox"/>	2002 <input type="checkbox"/>		Beige <input type="checkbox"/> Chestnut <input type="checkbox"/> Black <input type="checkbox"/> _____ <input type="checkbox"/> * Seasonal colors available

* To view seasonal colors, visit www.juzousa.com, then click on Product > Stocking > Fashion > Soft 2000