



Patient Last Name: _____ Patient First Name: _____
 Fitter Last Name: _____ Fitter First Name: _____
 Fitter Title: _____ (example: PT/OT/PTA)
 Date: _____

MEDI NON-CUSTOM FLAT KNIT SLEEVE MEASUREMENT FORM

Measure in centimeters

Left Right

Length Measurement from axilla to wrist (c-g) = _____



Measurement Points

Circumference **point 'g'** is taken around the upper arm at the axilla. For ease in determining this point, place a folded piece of paper up into the arm pit area. Point 'g' on the arm will be even with the top of the paper. **Measure with moderate tension.**

Circumference **point 'e'** is found in the middle of the elbow. **Measure without tension.**

Circumference **point 'c'** is the point of greatest compression and therefore the most important point. This point is found distally on the arm where the wrist meets the hand. One should find the smallest circumference at this point. **Measure without tension.**

Circumference 'g' = _____ **SILICONE BORDER ONLY**

Circumference 'e' = _____


Circumference 'c' = _____

Circumference 'a' = _____

Item	Class II only	Colors Available
MONDI ESPRIT	23 - 32 mmHg <input type="checkbox"/>	MEDI ESPRIT "CAMEL ONLY"

SIZE CHART

mondi esprit sleeves and gloves

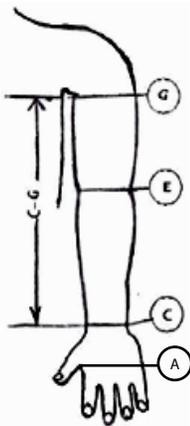

<i>circumference</i>	I	II	III	IV	V
a <i>HAND</i>	6 ³ / ₄ " - 7 ¹ / ₄ " 17.5-18.5 cm	7 ¹ / ₄ " - 7 ³ / ₄ " 18.5-20 cm	7 ³ / ₄ " - 8 ³ / ₄ " 20-22 cm	8 ³ / ₄ " - 9 ¹ / ₂ " 22-24 cm	9 ¹ / ₂ " - 10 ³ / ₄ " 24-27 cm
c <i>WRIST</i>	5 ³ / ₄ " - 6 ³ / ₄ " 15-17 cm	6 ³ / ₄ " - 7 ¹ / ₂ " 17-19 cm	7 ¹ / ₂ " - 8 ¹ / ₄ " 19-21 cm	8 ¹ / ₄ " - 9 ¹ / ₄ " 21-23 cm	9 ¹ / ₄ " - 10 ¹ / ₄ " 23-25 cm
e <i>ELBOW</i>	9 ³ / ₄ " - 11" 25-28 cm	10 ¹ / ₂ " - 11 ³ / ₄ " 27-30 cm	11 ¹ / ₂ " - 12 ¹ / ₂ " 29-32 cm	12 ¹ / ₄ " - 13 ¹ / ₄ " 31-34 cm	13" - 14 ¹ / ₄ " 33-36 cm
g <i>UPPER ARM</i>	11" - 12 ¹ / ₂ " 28-32 cm	11 ³ / ₄ " - 13 ¹ / ₂ " 30-34 cm	12 ¹ / ₂ " - 14 ¹ / ₄ " 32-37 cm	13 ³ / ₄ " - 15 ³ / ₄ " 35-40 cm	15" - 17 ¹ / ₄ " 38-44 cm
ARM LENGTH	d STANDARD: Less than 17 ¹ / ₄ " (44 cm) LONG: 17 ¹ / ₄ " (44 cm) and over				